

# Capillary Blood on Filter Paper for Determination of HbA<sub>1c</sub> by Ion Exchange Chromatography

JAN-OLOF JEPPSSON, MD, PHD  
PETER JERNTORP, MD, PHD  
LARS-OLOF ALMÉR, MD, PHD

RUTH PERSSON, RN  
GÖRAN EKBERG, MD  
GÖRAN SUNDKVIST, MD, PHD

**OBJECTIVE** — To facilitate HbA<sub>1c</sub> determination, we evaluated an HbA<sub>1c</sub> filter paper system enabling capillary blood sampling at home by the patients.

**RESEARCH DESIGN AND METHODS** — Capillary blood (two drops) was applied to a filter paper (HbA<sub>1c</sub> Via Post) and sent to the laboratory where a small disc was punched out on the filter paper. Hemoglobin was eluted from the disc in a buffer containing cysteine to eliminate the interfering glutathione adduct (HbA<sub>3</sub>) formed during storage. Analysis was performed by ion-exchange chromatography (Mono S, high-performance liquid chromatography), and the eluate was compared with hemolysate of venous blood from 41 patients. The stability of blood impregnated on filter paper was checked at different temperatures over different periods of time.

**RESULTS** — There was an excellent agreement ( $r = 0.99$ ) between HbA<sub>1c</sub> values from capillary blood on filter paper and HbA<sub>1c</sub> values from venous blood. HbA<sub>1c</sub> values were constant when stored on filter paper for 5–7 days at 20–21°C (room temperature) or at 4–6°C (refrigerator) for 10 days as well as at –70°C for several months after blood sampling. A new chromatographic-interfering hemoglobin fraction both from venous and capillary samples was identified as free  $\alpha$ -chain of hemoglobin.

**CONCLUSIONS** — The HbA<sub>1c</sub> filter paper system enables capillary blood sampling at home, eliminates the need of vein puncture in children and adults, and provides the diabetologist with an HbA<sub>1c</sub> value when the patient visits the clinic without a need for a previsit phlebotomy.

HbA<sub>1c</sub> has been used for almost two decades as a marker of long-term glycemic control in diabetic patients (1). The importance of HbA<sub>1c</sub> in this context has increased during the years, and HbA<sub>1c</sub> is now considered as the gold standard for metabolic control in diabetic patients; the improved glycemic control, as reflected by a near normalization of HbA<sub>1c</sub>, prevents the development of microvascular complications in type I diabetic patients (2,3). Indeed, without HbA<sub>1c</sub> determinations, it simply would not have been possible to demonstrate that improvement in metabolic control is associated with prevention of diabetic complications (4).

HbA<sub>1c</sub> is usually assayed from venous whole-blood samples. To enable analysis of capillary blood sampled on filter paper, a thiobarbituric acid colorimetric method and affinity chromatography have been introduced (5,6). When capillary blood was preserved and transported in liquid solutions, the correlation between venous and capillary HbA<sub>1c</sub> improved (7,8). Liquid specimens are difficult to handle for the patient; capillary blood samples on a filter paper are less complicated.

Assaying HbA<sub>1c</sub> from venous blood and mailing the HbA<sub>1c</sub> results directly to the diabetic patient improves glycemic control (9). We have now evaluated

whether capillary blood samples applied on a novel filter paper could be used in this context in combination with ion-exchange chromatography. In our HbA<sub>1c</sub> program, the patient takes a capillary blood sample at home and applies the blood to the filter paper. The filter paper is mailed to the laboratory, where the sample is eluted from the paper and HbA<sub>1c</sub> assayed. The aims of the current study were to optimize our conventional HbA<sub>1c</sub> method for the filter paper system, clarify essential technical steps in the procedure, relate capillary filter paper HbA<sub>1c</sub> values to HbA<sub>1c</sub> values obtained from venous whole blood, and to clarify whether HbA<sub>1c</sub> values are altered when dried spot blood samples on filter paper have been stored.

## RESEARCH DESIGN AND METHODS

### Blood samples

Venous blood was collected in EDTA-containing vacuum tubes (Becton Dickinson, Rutherford, NJ). Capillary blood (two drops) obtained from the fingertip was applied to a filter paper (Schleicher & Schuell, Keene, NH) attached to a request card: HbA<sub>1c</sub> Via Post (Boehringer Mannheim Scandinavia AB, Bromma, Sweden). The patient sampling was supervised by a nurse at the outpatient diabetes unit. The adult patients were selected at random at the same time they were venipunctured for ordinary HbA<sub>1c</sub> analysis. The dried capillary blood on the filter paper was stored in a small plastic bag and sent by regular mail (1–2 days in Sweden, ambient temperature 15–25°C) to the laboratory. If analysis was not performed directly on arrival to the laboratory, filter papers were maintained at 4°C (refrigerator) until analysis. To evaluate the effect of prolonged mailing, samples from 10 patients were sent around Sweden for 5 days at 15–25°C. The filter papers were eluted on the 6th day, and results were compared with venous blood from the 1st day when the capillary blood was applied on the filter paper.

From the Departments of Clinical Chemistry (J.-O.J.), Medicine (P.J., L.-O.A.), and Endocrinology (R.P., G.E., G.S.), University Hospital, Malmö, Sweden.

Address correspondence and reprint requests to Jan-Olof Jeppsson, MD, PhD, Department of Clinical Chemistry, University Hospital, S-20502, Malmö, Sweden.

Received for publication 24 January 1995 and accepted in revised form 14 September 1995.

HPLC, high-performance liquid chromatography.

### Sample preparation

**Venous blood.** Of the mixed whole blood, 17  $\mu\text{l}$  was diluted with 900  $\mu\text{l}$  sodium phosphate-citrate buffer also containing a detergent (0.02 and 0.05 mol/l, respectively, and 0.1% Triton X-100, pH 5.4). After incubation at 37°C for 30 min to enhance the elimination of pre-HbA<sub>1c</sub> (the labile fraction), the tubes were centrifuged, and 5  $\mu\text{l}$  was injected into the chromatographic system.

**Capillary blood on filter paper.** A small filter paper disc (3-mm diameter) was punched out from each paper on a spot with maximum amount of dried blood. Hemoglobin was eluted from the filter paper by adding 75  $\mu\text{l}$  of the sodium phosphate-citrate Triton X-100 buffer described above. After elution at 37°C during 30 min (to eliminate pre-HbA<sub>1c</sub>), an equal volume (75  $\mu\text{l}$ ) of a cysteine-containing buffer, pH 8.0 (Boehringer Mannheim Scandinavia AB), was added according to the manufacturer's instruction. After 45 min at room temperature, 5  $\mu\text{l}$  of the hemolysate was injected into the chromatographic system, and HbA<sub>1c</sub> analysis was performed as for venous blood.

### Equipment

The high-performance liquid chromatography (HPLC) system consisted of a pump (no. 2941, Pharmacia Biotechnology, Uppsala, Sweden), a Jasco 870 detector (Tokyo, Japan) equipped with a 415-nm filter, and a 10-mm flow cell together with a hydrogen lamp. The system contained an autoinjector for microtiter plates (BioRad, Richmond, CA) with a cooling system for 96 samples. A Shimadzu CR5A integrator (Tokyo, Japan) was used for calculating the peak areas using the valley-valley mode.

### Procedure

HbA<sub>1c</sub> was separated from other hemoglobin fractions by ion-exchange chromatography on a Mono S HR 5/5 column (Pharmacia, product no. 17-1040-01) by salt gradient elution at pH 5.7 (10). New Guard AX-300 (Applied Biosystems, Foster City, CA) was used as a precolumn. The flow rate was 2.0 ml/min, and the complete cycle time was 9 min. To prevent bacterial growth, sodium azide (0.02%, wt/vol) was added to all buffers, which were degassed under vacuum before use. If venous hemolysates and filter eluates were run on the same column, 100

$\mu\text{l}$  of the cysteine buffer solution was injected followed by a blank gradient before the patient samples.

### Stability

EDTA blood samples with different HbA<sub>1c</sub> concentrations from three patients were applied on multiple filter papers and stored at 20–21°C, 4–6°C, and –70°C for different time periods.

### Assay of imprecision

The between-day imprecision was determined by analyzing filter paper samples stored at –70°C during 6 months and by analyzing lyophilized blood samples (BioRad, Richmond, CA).

### Isolation of an interfering peak

One minor fraction between HbA<sub>1c</sub> and the major HbA<sub>0</sub> was pooled, dialyzed against water, and lyophilized. The  $\alpha$ - and  $\beta$ -chains were prepared, and their amino acid content and sequence were determined as previously reported (11).


## RESULTS

### Development of the method

Among the different fractions of hemoglobin appearing when dried blood is stored on filter paper, methemoglobin and a glutathione adduct of hemoglobin interfere with the chromatography pattern. Red cells contain a high amount of the tripeptide glutathione (glycine-cysteine-glutamic acid). The free thiol group of this peptide forms a disulfide bond with the  $\beta$ -chain of hemoglobin, and this compound increases during storage. The extracarboxyl group in glutamic acid contributes to an extra charged group in the hemoglobin molecule and forms a new hemoglobin fraction (HbA<sub>3</sub>) in the chromatographic profile (Fig. 1). The addition of a cysteine-containing buffer thereafter reduces the methemoglobin and eliminates the glutathione adduct. Hence, the chromatography profile after the addition of the cysteine-containing buffer was quite similar to that of native venous hemoglobin.

### Interpretation of chromatographic profile

The difference in net charge between the HbA<sub>1c</sub> fraction and other minor hemoglobin fractions may be small. Using the described salt gradient, a small fraction located to the right of HbA<sub>1c</sub> in the chromatography profile was noticed (X in Fig. 1).

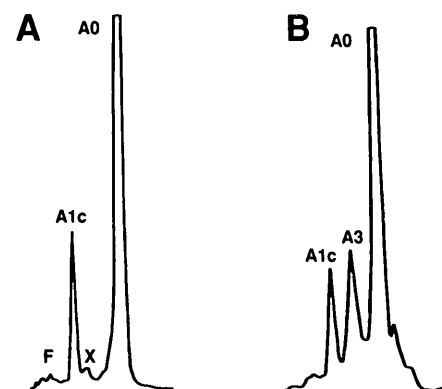


**Figure 1**—Ion-exchange chromatography of hemoglobin fractions from eluate of dried blood stored for 5 days on filter paper at room temperature: A) elution in the phosphate-citrate Triton X-100 buffer, pH 5.4, followed by addition of a cysteine-containing buffer, pH 8.0; B) results after elution in the phosphate-citrate Triton X-100 buffer. The fraction marked F represents the normal occurring fetal hemoglobin. X indicates free  $\alpha$ -chain that is separated from and not included in the HbA<sub>1c</sub> fraction. HbA<sub>3</sub> is the glutathione adduct formed during storage both in filter paper and aged venous EDTA blood.

Amino acid analysis and sequence determination of this fraction (in nondiabetic and diabetic subjects) demonstrated an identity with purified  $\alpha$ -chain, indicating that this interfering fraction contained dimers or tetramers of the  $\alpha$ -chain of hemoglobin. The fraction was visible in chromatography patterns both from fresh venous blood and filter paper eluate.

### Comparison between venous blood and capillary blood

HbA<sub>1c</sub> from venous blood and capillary blood sampled on the same occasion from 41 patients were compared using a similar chromatography system for analysis of HbA<sub>1c</sub> in venous and capillary blood. Figure 2 shows the close correlation between HbA<sub>1c</sub> values obtained from venous blood analyzed on day 1 (when the blood sample was taken) and HbA<sub>1c</sub> values obtained from dried capillary blood samples (taken simultaneously with the venous blood sample) but stored on filter paper at room temperature in darkness and kept in plastic bags for 2 days before the assay. The correlation was as close for low values of HbA<sub>1c</sub> as for high values. The mean value and range for filter paper and venous blood was 7.49% (4.37–12.28) and 7.33% (4.13–12.21), respectively. Filter paper gave an average of 0.16% higher values than did venous blood in this study.

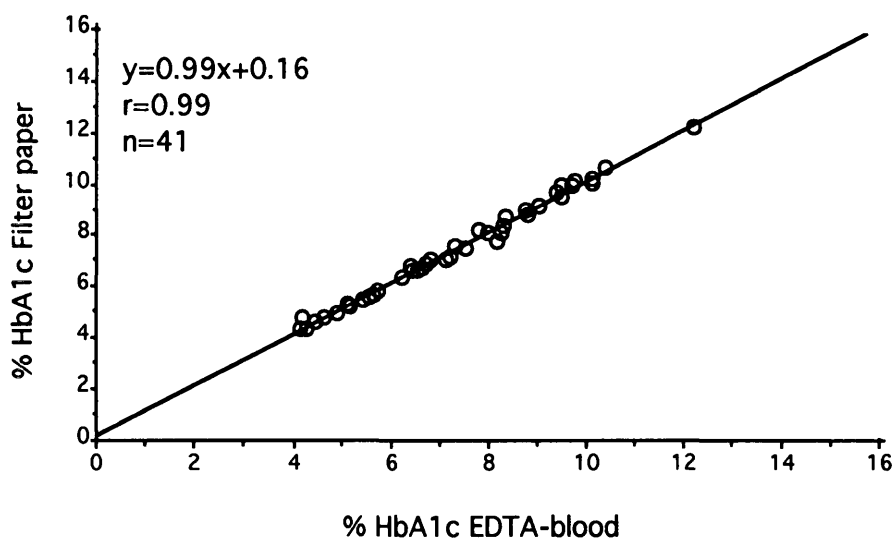

**Figure 2**—Comparison of HbA<sub>1c</sub> obtained from filter paper (capillary samples) and EDTA blood (venous whole blood) on a Mono S column for HbA<sub>1c</sub> on a standard HPLC system.


**Storage and stability of capillary samples**

Samples from individual patients showed no significant differences in HbA<sub>1c</sub> results when stored for 5 days at room temperature (in darkness) or 4–6°C for 10 days and for several months at –70°C after blood sampling (Fig. 3).

Mailing for 5 days during a weather period with a daily outside temperature of 15–25°C did not produce any change in the chromatographic profile. The between-day imprecision (coefficient of variation) of venous and capillary blood samples within the reference range was 1.35 and 1.54%, respectively.

**CONCLUSIONS** — In this study, we demonstrate that capillary-sampled blood collected on filter paper shows HbA<sub>1c</sub> values corresponding to those obtained by conventional venous whole-blood samples. Moreover, storage of capillary samples on filter paper at 20–21°C for 5 days did not change the HbA<sub>1c</sub> results. The tested HbA<sub>1c</sub> system therefore seems appropriate for use in diabetic patients. This system will allow the diabetic patient to take a capillary blood sample at home and then send the sample by mail to the laboratory for a reliable HbA<sub>1c</sub> determination. Two years of experience with frequent comparisons between venous and capillary blood samples show that the introduced technique works in a home setting if the patient is orally informed by a nurse at the first sampling. Every request form also contains a clear explanation

and instruction to the patient. HbA<sub>1c</sub> from both venous blood and filter paper eluates was quantitated after a minor modification (addition of sodium azide)

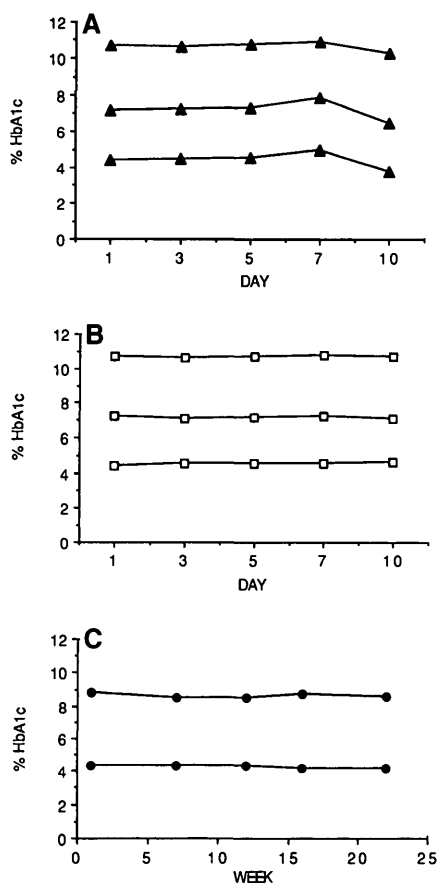

**Figure 3**—Storage of filter paper samples in A) room temperature (20–21°C), B) refrigerator (4–6°C), and C) freezer (–70°C).

of our original ion-exchange chromatography method (12). The modified system had sufficient resolution to give accurate results for the filter paper eluate.

The 97.5 percentile value of HbA<sub>1c</sub> from healthy individuals in this study was 4.9% and in good agreement with Bisse and Wieland (13) but considerably lower than that reported by Little et al. (14). It is well known, however, that the same blood specimen tested in different HbA<sub>1c</sub> systems, such as ion-exchange chromatography, affinity chromatography, electrophoresis, or mono- and polyclonal antibody techniques, produces different results. To avoid confusion, an international calibrator for HbA<sub>1c</sub> has to be introduced (15). Nevertheless, capillary blood samples on filter paper have the potential to be used in most HbA<sub>1c</sub> assays, including immunoassay (16), and not only HPLC, as in this study.

A filter paper-based HbA<sub>1c</sub> system allows HbA<sub>1c</sub> results to be available when patients visit their doctor. At our hospital, >2,000 capillary samples have now been handled in this way, eliminating previsit venipunctures. Besides reducing the need for the patient to come to the clinic for phlebotomy, capillary HbA<sub>1c</sub> also diminishes the burden for the diabetic team; the need for venipunctures by diabetic nurses or laboratory technicians are substantially reduced. The gain in the diabetic child is obvious; fear of phlebotomy that might interfere with the goal in diabetic care is avoided. Intensive attention improves glycemic control (17), and this is promoted by mail delivery of HbA<sub>1c</sub> results directly to the diabetic patient at home (9). This is easily accomplished by the capillary blood HbA<sub>1c</sub> filter paper system. To further facilitate for the diabetic patient, we provide our HbA<sub>1c</sub> results directly in a graphic form from the laboratory computer presenting the current and previous HbA<sub>1c</sub> test results (absolute values together with an HbA<sub>1c</sub> graph). The form is mailed to the patient simultaneously as the result is sent to the patient's diabetologist. Another advantage with this HbA<sub>1c</sub> system is the quality. Submitting samples by mail will allow laboratories to conduct adequate HbA<sub>1c</sub> assays with the same accuracy and precision from both diabetic wards and outpatient units.

Certainly, after the Diabetes Control and Complications Trial (2), HbA<sub>1c</sub> will remain a focus for many years. A

normal or slightly deviated HbA<sub>1c</sub> value without the occurrence of major adverse hypoglycemic occasions has to be the characteristic of diabetes care in the future. To promote this, an HbA<sub>1c</sub> system where the patient takes a capillary blood sample and mails it to a laboratory that uses a valid HbA<sub>1c</sub> assay must take place. The result may not only be provided to the diabetologist but also to the patient on an appropriate and didactic form.

**Acknowledgments**—The filter papers and reagents were provided by Boehringer Mannheim Scandinavia AB, Bromma, Sweden.

The skillful assistance of Benny Larsson and Anna Arnetorp is gratefully acknowledged.

### References

1. Goldstein DE: Is glycosylated hemoglobin useful? *N Engl J Med* 310:284–285, 1984
2. The Diabetes Control and Complications Trial Research Group: The effect of intensive treatment of diabetes on the development and progression of long-term complications in insulin-dependent diabetes mellitus. *N Engl J Med* 329:977–986, 1993
3. Reichard P, Nilsson BY, Rosenqvist U: The effect of long-term intensified insulin treatment on the development of microvascular complications of diabetes mellitus. *N Engl J Med* 329:304–309, 1993
4. Raskin P: The DCCT: finally. *J Diabetes Complications* 7:214–215, 1993
5. Goldstein DE, Wiedmeyer HM, England JD, Little RR, Parker KM: Glycosylated protein in whole blood spotted on filter paper. *Clin Chem* 28:386–387, 1982
6. Little RR, McKenzie EM, Wiedmeyer HM, England JD, Goldstein DE: Collection of blood on filter paper for measurement of glycosylated hemoglobin by affinity chromatography. *Clin Chem* 32:869–871, 1986
7. Voss EM, Cembrowski GS, Clasen BL, Spencer ML, Ainslie MB, Haig B: Evaluation of capillary collection systems for HbA<sub>1c</sub> specimens. *Diabetes Care* 15:700–701, 1992
8. Eckerbom S, Bergquist Y: Improved sample collection technique for capillary blood on filter paper for determination of glycosylated hemoglobin. *Ann Clin Biochem* 26:148–150, 1989
9. Jerntorp P, Jeppsson JO, Sundkvist G: Hemoglobin A<sub>1c</sub> self-recording in the management of diabetes mellitus. *Acta Med Scand* 223:359–363, 1988
10. Jeppsson JO, Jerntorp P, Sundkvist G, Englund H, Nylund V: Measurement of hemoglobin A<sub>1c</sub> by a new liquid chromatography assay. *Clin Chem* 32:1867–1872, 1986
11. Jeppsson JO, Källman I, Lindgren G, Fågerström L: Hb-linköping: a new hemoglobin mutant characterised by reversed-phase high-performance liquid chromatography. *J Chrom* 297:31–36, 1984
12. Eckerbom S, Bergquist Y, Jeppsson JO: Improved method for analysis of glycosylated haemoglobin by ion exchange chromatography. *Ann Clin Biochem* 31:161–167, 1994
13. Bisse E, Wieland H: High-performance liquid chromatography separation of human haemoglobins. *J Chrom* 434:95–110, 1988
14. Little RR, Wiedmeyer H-M, England JD, Wilke AL, Rohlfing CL, Wians FH Jr, Jacobson JM, Zellmer V, Goldstein DE: Inter laboratory standardisation of measurements glycohemoglobins. *Clin Chem* 38:2472–2478, 1992
15. Bruns DE: Standardisation calibration, and the care of diabetic patients. *Clin Chem* 38:2363–2364, 1992
16. Jeppsson JO: Determination of HbA<sub>1c</sub> by the Tina-quant HbA<sub>1c</sub> immunoassay using dried capillary blood on filter paper. *Klin Lab* 39:1080–1082, 1993
17. Worth RC, Home PD, Johnston DG, Anderson J, Ashworth L, Barrin JM, Appleton D, Binder C, Alberti KG: Intensive attention improves glycemic control in insulin-dependent diabetes without further advantage from home blood glucose monitoring: results of a controlled trial. *Br Med J* 285:1233–1240, 1982